

The NEXT-DEMO++ detector

F. Monrabal,^a I. Liubarsky,^a S. Carcel,^a J.J. Gómez-Cadenas,^a

^a*Instituto de Física Corpuscular (IFIC), CSIC & Universidad de Valencia
46980 Valencia, Spain*

E-mail: francesc.monrabal@ific.uv.es, Igor.Liubarsky@ific.uv.es,
carcel@ific.uv.es, gomez@mail.cern.ch

Contents

1	Introduction	2
2	Pressure vessel	3
3	Gas system	4
3.1	Cryo-recovery protocol	7
3.2	Gas Mixtures	7
3.3	Technical and safety specifications	8
4	Field Cage	8
4.1	Drift region	9
4.2	Buffer region	9
4.3	High voltage feedthroughs	9
4.4	Cathode grid	10
4.5	Electroluminescent region	11
5	Electrical Power	12
6	Experimental area distribution	14
7	Safety issues	16
7.1	Pressure Risks	16
7.1.1	Pressure Vessel	16
7.1.2	Gas system	16
7.1.3	Recovery Bottles	16
7.2	Oxygen displacement	16
7.2.1	Safety systems implemented	17
7.3	High Voltage safety	17
7.4	Fire	17
7.5	Cryogenics	17
8	Risk assesment for gas	18
8.1	Characteristics of the area	18
8.2	Related Gas systems	21
8.3	Gas release risk assesment	21
8.3.1	Relevant Facts	21
8.3.2	Risk assesment	21
8.4	Conclusions	21
8.4.1	Oxygen detection hazard	21

1 Introduction

NEXT-DEMO is a high-pressure xenon TPC contained within a cylindrical stainless-steel pressure vessel of diameter 30 cm and length 60 cm which was designed to withstand up to 10 bar. The TPC itself is defined by three metallic wire grids — called *cathode*, *gate* and *anode* — which define the two active regions: the 30-cm long *drift region*, between cathode and gate with a drift field up to 800 V cm^{-1} ; and the 5 mm long *EL region*, between gate and anode. The electric field is created by supplying a large negative voltage to the cathode, then degrading it using a series of metallic rings of 30 cm diameter spaced 5 mm and connected via $0.5 \text{ G}\Omega$ resistors (shown in figure 1). The gate is at negative voltage so that a moderate electric field — $[1.0, 4.0] \text{ kV cm}^{-1} \text{ bar}^{-1}$ — is created between the gate and the anode, which is at ground. A set of six panels made of PTFE (Teflon) coated with tetraphenyl-butadiene (TPB) are mounted inside the electric-field cage forming a *light tube* of hexagonal cross section with an apothem length of 8 cm.

The NEXT-DEMO detector has been operating continuously at IFIC (Figure 2 shows the NEXT-DEMO detector in the experimental are at IFIC)without major problems

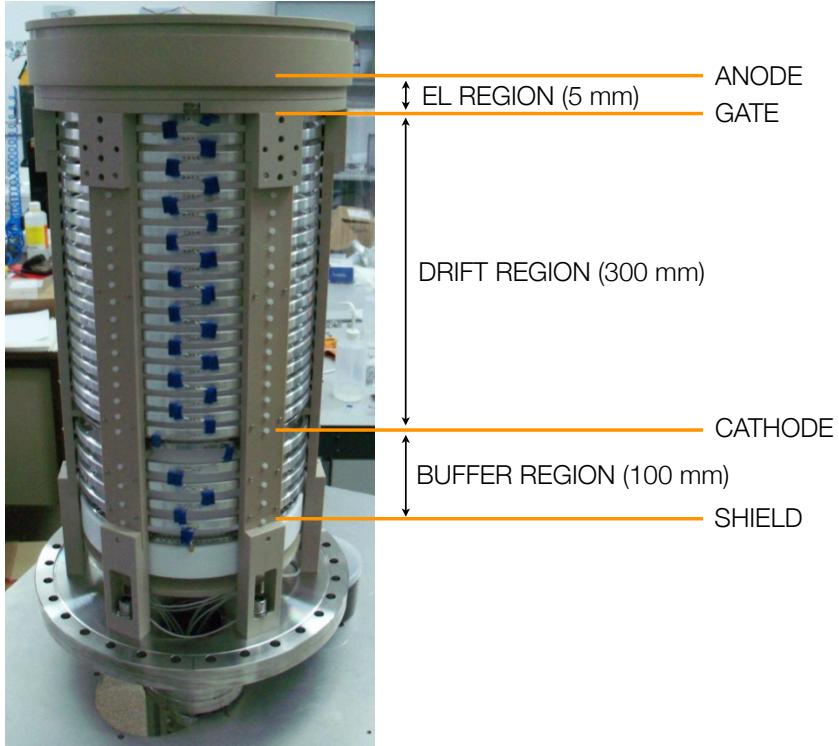


Figure 1. External view of the time projection chamber mounted on one end-cap. The approximate positions of the different regions of the TPC are indicated.



Figure 2. Picture of the NEXT-DEMO detector in the experimental area at IFIC with all the different systems needed for its operation: Gas system, High voltage modules and electronics.

during the last three years (2011-2014). Its operation and results has been useful to fully demonstrate the capabilities of the NEXT technology for the search of the neutrinoless double beta decay process.

The next step for the NEXT-DEMO detector is to keep helping as a demonstrator for new ideas that can be implemented in larger detectors. In that sense, the possibility to use a magnetic field in parallel to the drift field of the TPC could give us an improvement in the background rejection factor of almost an order of magnitude with no significant lose of efficiency. This affirmation comes from our recent Monte-Carlo studies and it needs to be tested experimentally with a realistic set-up. Here the NEXT-DEMO detector seems to be the clear candidate.

The project that we propose to CERN consists in placing the NEXT-DEMO detector inside the HARP magnet at CERN and operate it in order to confirm the Monte-Carlo results.

2 Pressure vessel

The pressure vessel of NEXT-DEMO, shown in figure 3, is a stainless-steel (grade 304L) cylindrical shell, 3 mm thick, 30 cm diameter and 60 cm length, welded to CF flanges on both ends. The two end-caps are 3-cm thick plates with standard CF knife-edge flanges. Flat copper gaskets are used as sealing. The vessel was certified to 10 bar operational pressure. It was designed at IFIC and built by Trinos Vacuum Systems, a local manufacturer. Additional improvements — including the support structure and a rail system to open and move the end-caps — have been made using the mechanical workshop at IFIC.

The side of the chamber includes 8 CF40 half-nipples. One set of 4 is located in the horizontal plane while the other is displaced towards the underside with respect to the

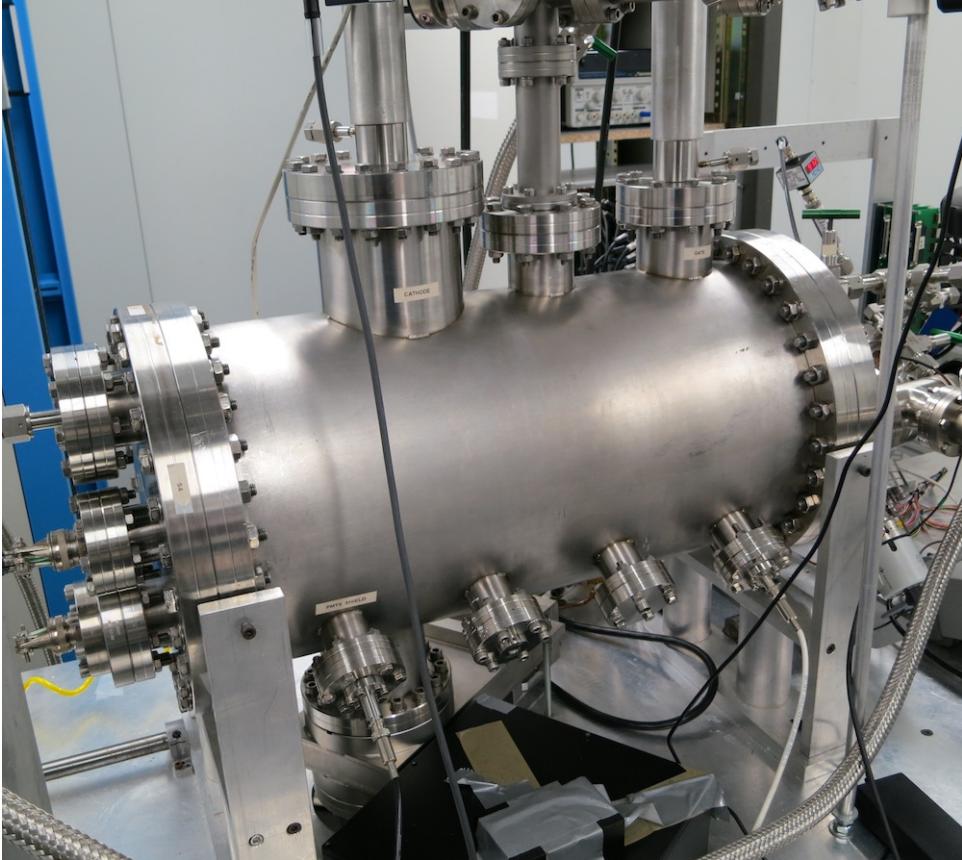


Figure 3. The pressure vessel of NEXT-DEMO.

first set by 60° . These contain radioactive source ports used for calibration of the TPC. The ports are made by welding a 0.5 mm thick blank SS plate onto a 12 mm OD pipe on a CF40 liquid feedthrough. On top of the vessel and along the vertical plane there are three additional half-nipples (CF130, CF67 and CF80) used for high-voltage input and connection to a mass spectrometer (through a leak valve). On the opposite side, at the bottom, a CF100 port connects the pressure vessel to the vacuum pumping system. A guillotine valve closes this connection when the vessel is under pressure. The end-caps include several CF ports for the connections to the gas recirculation loop and for the feedthroughs (power and signal) of the PMT planes.

3 Gas system

The functions of the gas system of NEXT-DEMO are the evacuation of the detector, its pressurization and depressurization with xenon (and argon), and the recirculation of the gas through purification filters. A schematic of the system is shown in figure 4.

The standard procedure during normal operation of the detector starts with the evacuation of the vessel to vacuum levels around 10^{-5} mbar. The detector is then filled with xenon gas to pressures up to 10 bar. The xenon can be cryogenically recovered to a

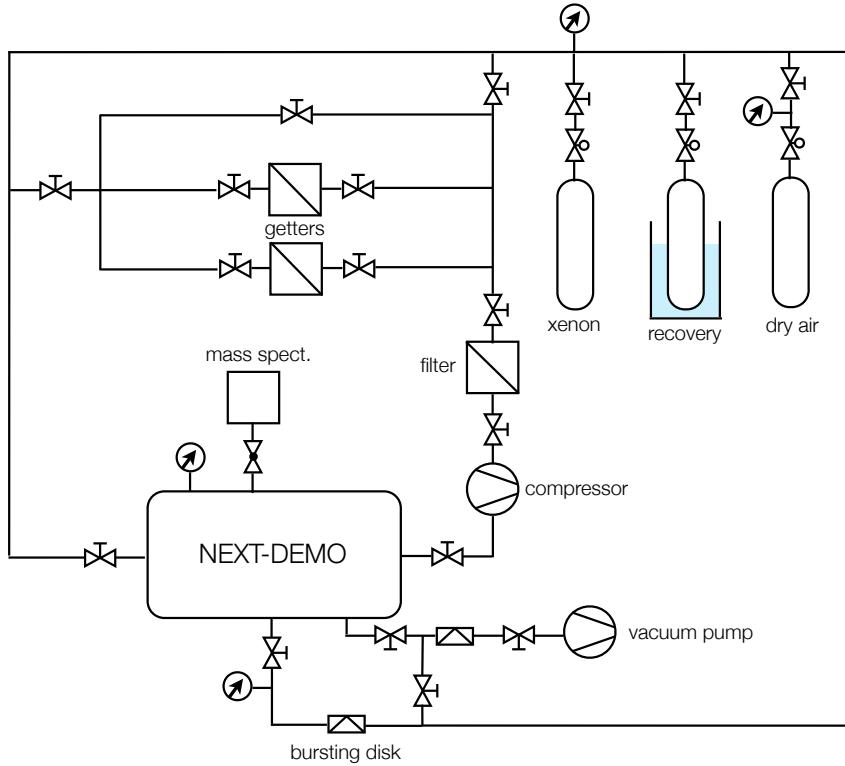


Figure 4. Simplified schematic of the gas system of NEXT-DEMO.

stainless-steel bottle (Fig. 5) connected to the gas system by simply immersing this in a dewar filled with liquid nitrogen. The gas flows inside the bottles and freezes there.

The vacuum pumping system consists of a roughing pump (Edwards XDS5 scroll vacuum pump) and a turbo molecular pump (Pfeiffer HiPace 300). Vacuum pressures better than 10^{-7} mbar have been obtained after pumping out the detector for several days. The recirculation loop is powered by an oil-less, single-diaphragm compressor (KNF PJ24999-2400) with a nominal flow of 100 standard liters per minute. This translates to an approximate flow of 10 liters per minute at 10 bar, thus recirculating the full volume of NEXT-DEMO (~ 45 L) in about 5 minutes. The gas system is equipped with both room-temperature (SAES MC50) and heated *getters* (SAES PS4-MT15) that remove electronegative impurities (O_2 , H_2O , etc.) from the xenon. All the gas piping, save for the inlet gas hoses and getter fittings, are 1/2 inch diameter with VCR fittings. A set of pressure relief valves (with different settings for the various parts of the system) and a burst disk in the vacuum system protect the equipment and personnel from overpressure hazards.

The operation of the gas system (Fig. 6) has been, in general, very stable. The detector has run without interruption for long periods with no leaks and continuous purification of the gas. However, one major leak occurred when the diaphragm of the recirculation pump broke, causing the loss of the xenon volume contained in the chamber. This led to the

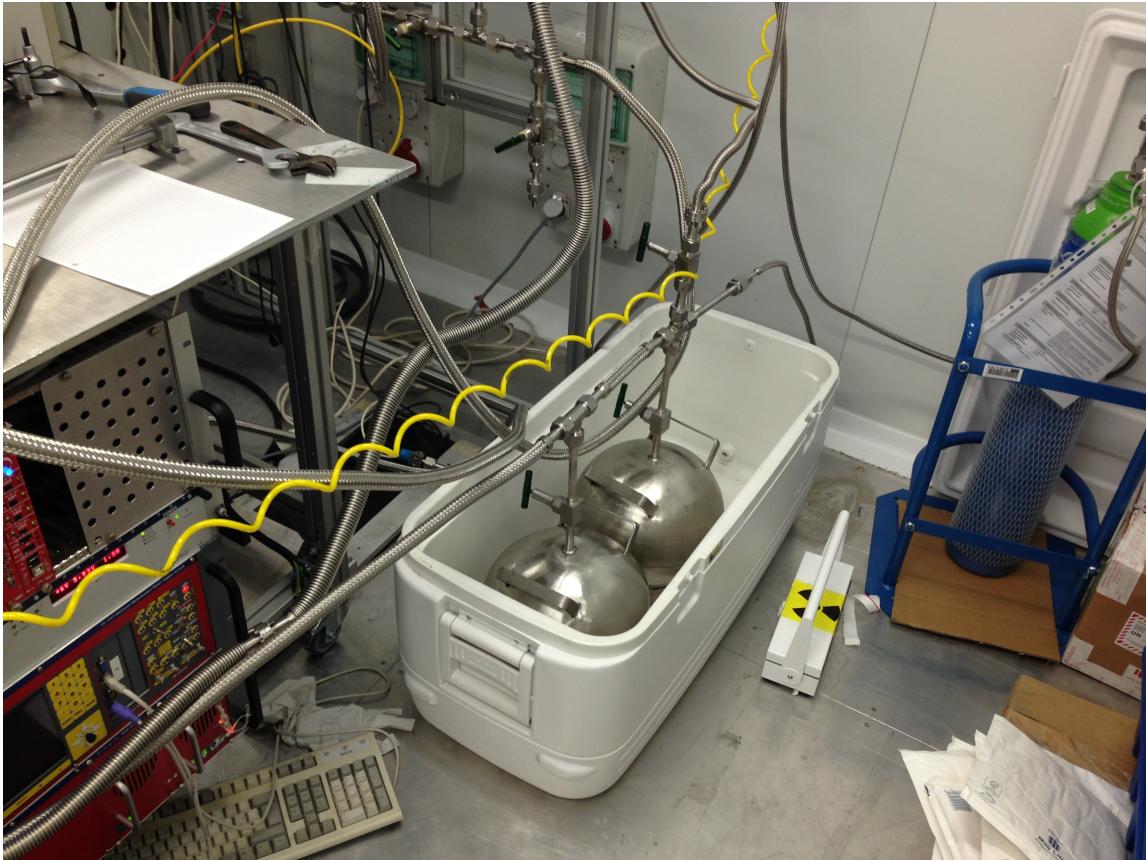


Figure 5. Stainless steel bottles used for cryogenic recovery of the Xenon

installation of an emergency mechanism that, in the event of pressure drop, automatically closes those valves connecting the pump to the rest of the gas system. Since installation, only one major failure has taken place which the emergency system isolated without loss of gas. Micro-leaks, on the level of 0.005 bar per day, due to bad connections in the gas system have also been detected making it necessary to introduce additional xenon to maintain pressure. The micro-leaks were found and properly repaired allowing for a stable detector operation in the last year.

Several improvements to the initial design of the gas system have been made thanks to the initial data runs. These include the recognition of the importance of reliability of the main pump as well as the decision to use hot getters as the main gas purification stage due to their negligible emission of radon compared to that of room-temperature getters.



Figure 6. Picture of the Gas system used in the NEXT-DEMO detector.

3.1 Cryo-recovery protocol

When the gas system needs to be stopped for an intervention in the detector the gas Xenon needs to be recovered and saved. The way of recovering the Xenon from the vessel and the gas system is with cryogenics.

NEXT-DEMO gas system has two spherical bottles specially designed for this use (see attachments) that are connected to the vessel and to the gas system. The bottles are placed in a thermally insulated container that is filled with liquid nitrogen until it reaches (approx.) half of the height of the bottles. Once the nitrogen stops boiling we open the corresponding valve to the part of the system that we want to recover, vessel or gas system. It is also possible to have both open at the same time if needed. Then the Xenon flows into the recovery bottles where it freezes.

Once the pressure in the vessel/gas system is stable at 0.03 absolute bar the bottles are full and we close them. Then we need to wait for the liquid nitrogen to evaporate.

3.2 Gas Mixtures

For the upgrade of NEXT-DEMO++ we will need to operate the detector with a mixture of Xenon with other gas. The gas chosen as a first option is CO₂. This gas give us all the advantages of the gas mixtures (improvement in drift velocity, reduction of diffusion,...)

at concentrations smaller than 1% . In order to control the concentration of CO₂ the gas system will need a small upgrade consisting in two flow meters and a small bottle of CO₂.

3.3 Technical and safety specifications

The gas system is made up of standard components. The Getters and the compressor carry CE marks.

Swagelock components do not carry a CE mark. However, all Swagelock components used comply with the requirement of Article 3, Paragraph 3 of the PED Directive (97/23/EC) and in accordance with that section do not carry a CE mark. Notwithstanding, we can included all the certificates.

The complete list of parts with the associated part numbers is given in attachment "VALCI-MONTI-195". We have also selected a list of major components that is given below, in which we have omitted explicit listing of the standard components such as elbows and similar coupling as they are too numerous and being all manufactured by Swagelock comply with the above mentioned PED directive.

- 28 × 1/2 VCR valves SS/8BG-VCR
- 3 × 1/4 VCR valves 6LVV-DPHVR4-P1
- 2 × 1/4 VCR valves SS-4BC-VCR
- 3 × regulators KPR1JRFX27A20000
- 2 × Cold getters SAES Pure Gas MC450-902
- 1 × Hot getter SASE Pure Gas PS4-MT15-R2
- Compressor: 1 × KNF PJ24999-2400

Non standard components are: The Vessel and the Stainless steel recovery bottles which were specifically designed and tested for NEXT. The calculations and test document are included in the attachment document to this one.

4 Field Cage

The NEXT-DEMO field cage is designed to provide an homogeneous and uniform electric field inside the active volume of the NEXT-DEMO detector. To achieve this different parts have to be designed, constructed and assembled. The different components of the the field cage project are:

- Drift region.
- Buffer.
- High voltage feedthroughs.
- Cathode grid.

- Electroluminescent region.

The main safety issues are related with the high voltage that we need to apply to the detector. In that sense, the entire detector and all of its metallic components are to be grounded for safety. The grounding points will provide path of least resistance to any stray transients resulting from HV discharges as is common for operation of HV equipment. Moreover, the currents will be small ($\mathcal{O}(\mu A)$) due to the large value of the resistors involved.

4.1 Drift region

The drift regions needs to create a moderate (300-600V/cm) but very homogeneous electric field. This is done using copper rings connected $10G\Omega$ resistors. Figure 7 shows an example of the rings used for the NEW detector, the concept used in DEMO will be the same.

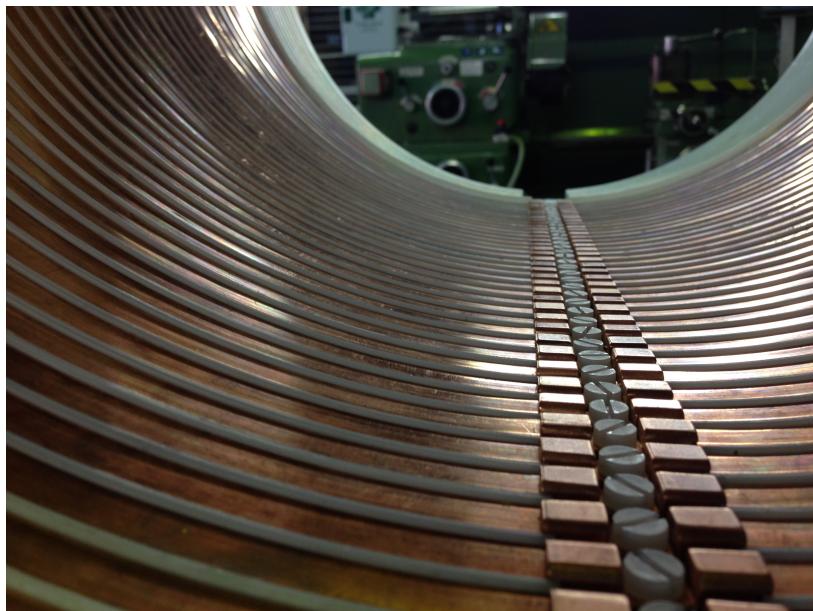


Figure 7. Detail of the copper rings in the drift region.

4.2 Buffer region

In order to protect the PMTs' windows from the cathode voltage safely a certain distance is needed. This space is what we call buffer region. In the current design the buffer region consists in a structure made completely with High Density Polyethilene (HDPE). The electric field in this region is much higher than in the drift region, as much as 5kV/cm, and it can result in sparks in different parts of the detector. For that reason the buffer has been extensively tested and has demonstrated to be very robust and stable.

4.3 High voltage feedthroughs

In order to create the different field regions high voltage has to be applied in different parts of the detector. Due to the impossibility of generating the high voltage (HHV) inside the

detector (due to radioactivity and space problems), the HHV has to be produced outside and then passed inside using feedthroughs that passes the voltage but stops gas leaks.

When discussing about feedthroughs we should differentiate three parts: The high voltage modules, the cable and the feedthrough itself.

Figure 8) shows the design and the results of the electric field simulations.

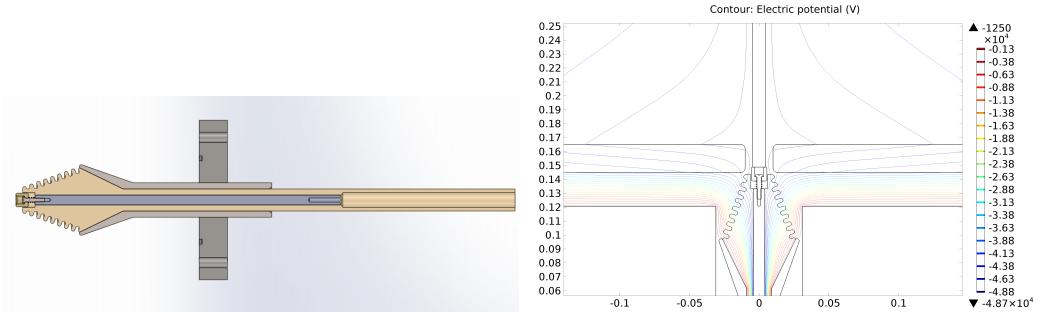


Figure 8. Design (left) and Comsol simulation (right) of the HVFT to be used in NEXT-DEMO.

A key part on that design was the cryo-fitting components to create gas tight seals as all different parts of the HVFT are sealed using this technique.

In the other hand, a few considerations need to be taken into account when operating the detector:

- The HHV modules need to be grounded using a proper connexion to the rack frame. This connexion should be checked periodically to ensure the quality of the connexion.
- The cable needs to be supported with a flexible/movable structure at the top of the lead castle that will allow to hold the cable when opening/closing the castle without damaging the cable.

4.4 Cathode grid

The cathode grid consists of a stainless steel frame with wires to fix the potential. Figure 9 shows a detail of the cathode frame with the grooves for fixing the wires and the bronze tensioners.

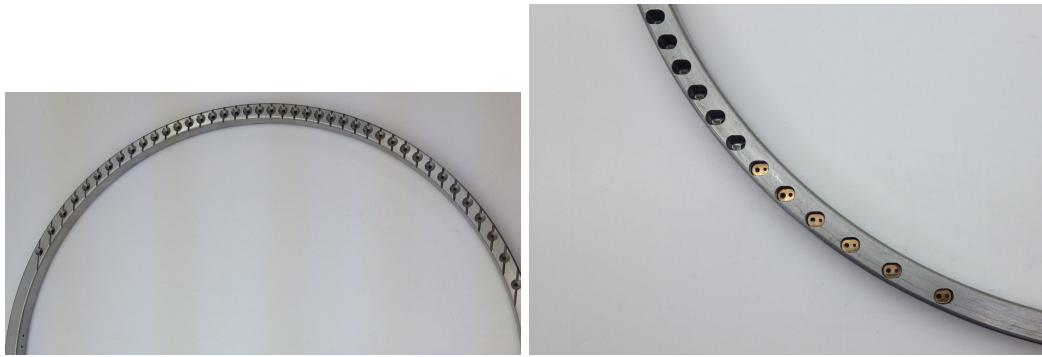


Figure 9. Cathode frame (left) with a detail of the groves for fixing the wires. The bronze tensioners (left) will be used to strength avery wire at the right tension.

4.5 Electroluminescent region

The electroluminescent region is the amplification region of the detector and it is designed to work with the highest electric field. It consist of a steel mesh grid and a fused silica plate both supported by a HDPE frame.

The main parts of the Electroluminescent region are the mesh and the anode plate. The anode consists in a plate of fused silica coated with ITO in one side and TPB in the other.

Figure 10 shows details of the different parts of the electroluminescent region.



Figure 10. Top pictures show the different frames used for assembly and tension the gate mesh. Bottom is a picture of the fused silica plate before coating with ITO

5 Electrical Power

The major part of the electrical power needed for the operation of NEXT-DEMO detector is related with the recirculation pump. The other components of the gas system will only require a fraction of that power. Table 1 reflects the power consumption by the components of the gas system.

About the power needed for the operation of the HV and the electronics, table 2 reflects the power consumption of the different systems as the total power needed.

In order to operate in a safer mode the final power request will be about 40% more than the nominal power consumption. In that case the request will be 15.8kW.

Recirculation Pump PC	3kW
Roughing pump	250 W
Hot getter	660 W
Total	3.91 kW

Table 1. Power consumption for NEXT-DEMO++ gas system

DAQ PC	5kW
SiPM FE Boards	1.5kW
FECs	500 W
ATCA	300W
SiPMs	30W
Total	7.33 kW

Table 2. Power consumption for NEXT-DEMO++ electronics

The power supply sources used for the electronics are

6 Experimental area distribution

The experimental area provided for the operation of the NEXT-DEMO detector is $25m^2$, thus we have assumed a square of 5x5m. In this area we have made the best possible distribution that allow us to have everything in the room and, at the same time, give us enough clearance to move inside the room safely.

Figure 11 shows the proposed distribution of NEXT-DEMO, the magnet, and all the different parts needed for its operation. The details of this distribution need to be discussed with CERN.

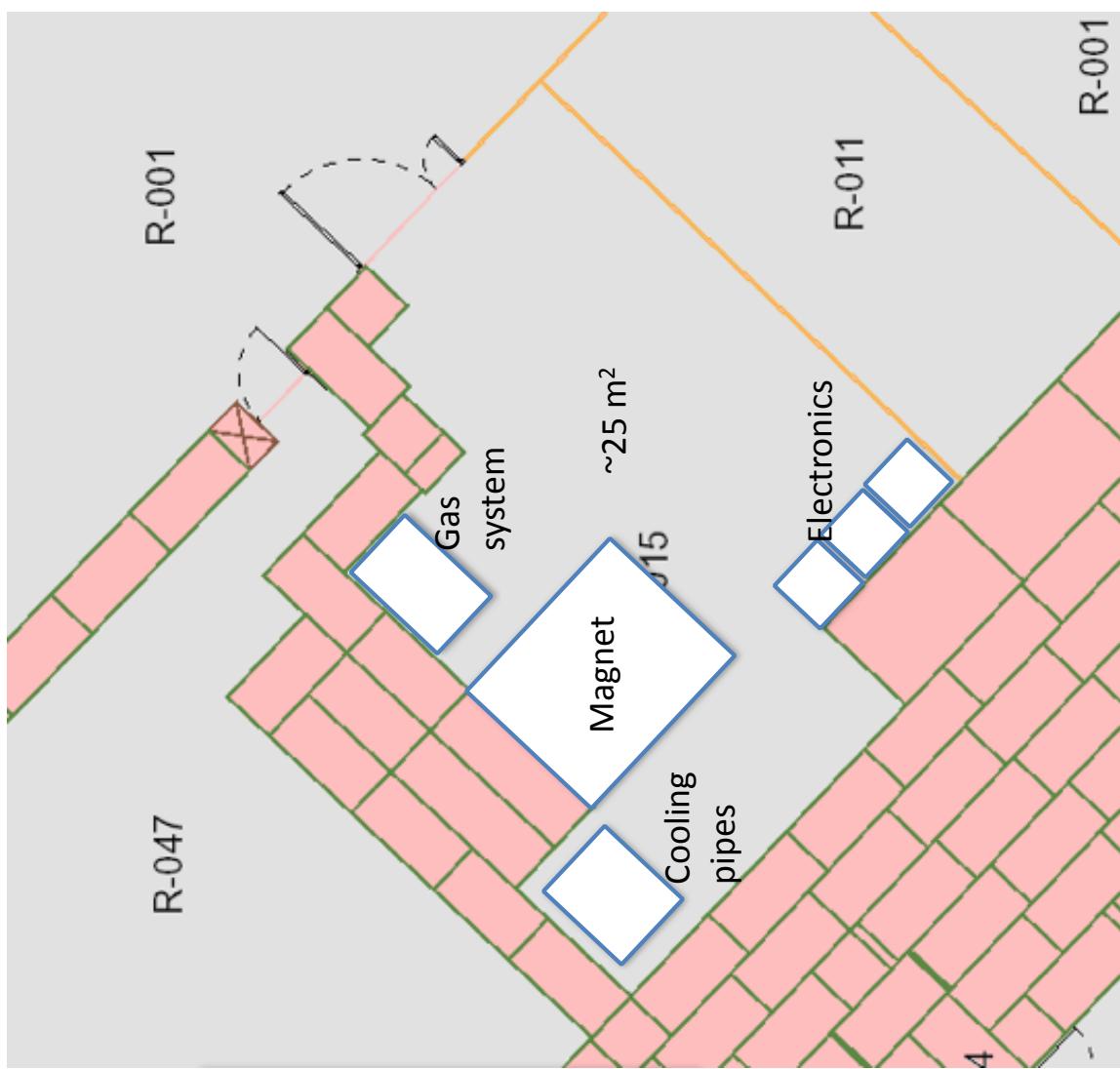


Figure 11. Proposed distribution of the NEXT-DEMO, magnet and all different necessary systems for the operation of the detector in the experimental area

7 Safety issues

In this section we review the main issues related with safety when operating and working with the NEXT-DEMO detector evaluating the risk level for all of them. We do not include the issues related with the operation of the magnet.

7.1 Pressure Risks

The components that are related with pressure risks are the pressure vessel and the gas system. Both have been operating for about three years with no major problems.

7.1.1 Pressure Vessel

The engineering calculations and the hydrostatic test results on the vessel are attached to this document. These documents show that the operation of the NEXT-DEMO pressure vessel is safe and it doesn't represent a risk for its operation.

7.1.2 Gas system

The gas system is composed of standard swagelock components, all of them with the european certification for safety (CE mark). The list and scheme of the different components is shown in the correspondind attached documents. The gas system does not represent a significant risk for the operation of the NEXT-DEMO detector.

7.1.3 Recovery Bottles

The bottles designed for the cryogenic recovery of the Xenon from the main volume have been tested for pressure and they are safe for operation. The document of such test are included in the attahments.

7.2 Oxygen displacement

When operating with any gas we need to understand the effect that a total failure of the system will have on the amount of oxygen in the experimental area. The amount of Xenon used in the NEXT-DEMO detector is around 1kg, for our safety estimations we will assume 1.5kg so we are in a very pessimistic scenario. In case of a big leak in gas system all the Xenon will be liberated into the experimental area. Operating at 10 bar and assuming a leak with a size of a 1/4" inch hole, the flow will be of the order of $0.5 \text{ m}^3/\text{min}$. The volume of 1.5kg of Xenon at 1 bar represents a volume of 8.84 cubic meters, it will take about 15 miutes to release the whole volume to the atmosphere. As Xenon is much denser than air it will stay at the bottom side. The height of the volume that Xenon will occupy is defined by the total amount of Xenon and the surface of the experimental area. As the surface of the experimental area is 25 square meters, the total height occupied by the Xenon will be 0.35 meters assuming the experimental area is sealed. That should not represent any danger for working in the area but oxygen monitors, either individual or general, are recomendend to reduce this risk. Also, the slow control system switch off the pump and locks and send an e-mail giving information of a pressure drop. A light alarm can be installed in the experimental area and connected to the slow control if necessary.

7.2.1 Safety systems implemented

A slow control system has been developed in order to monitor the pressure in different parts of the gas system and pressure vessel. The slow control can start/stop the re-circulation pump and also block the system in case a leak is detected so the gas remains inside the system. The slow control represents an extra help for improving the safety.

Moreover, the gas system has two bursting disks for a controled break in case of over-pressure. One in the vacuum side for protecting the vacuum instruments (turbo molecular pump, RGA,...) and one in the pressure side that breaks at 13 bar to prevent any high over-pressure of the system.

7.3 High Voltage safety

The high voltage to create the electric field in TPC is generated by two High voltages modules (*FUG HCP 140-100000* and *FUG HCP 0-35000*). Those modules have a fine tunning potentiometer for both voltage and current, allowing us to limit the total power and current produced by the module. Moreover, they can be controled remotely using an IP connexion. During normal operation the modules are always controled using a special Slow Control software developed by the IFIC grup using LabView[©] libraries. That prevents any errors in the manipulation of the modules and also helps with fast response in case of any failure. The documentation related with the modules can be found in the folder "HHV Power Supplies" Moreover, there is no direct acces to the any point that is at high voltage, all of them will be either inside the detector or inside the modules. That shows that the high voltage operation does not represent a risk for the operation of the detector.

7.4 Fire

As the gases used in the operation of NEXT-DEMO are not flamable the only system that is affected by the risk of fire is the electronics. The total power is not very high (see section 5) and the electronics has been operating for a long time without any problem, we consider that the risk of having fire in the experimental area is minimum. In case that a extinguisher system is needed we can design a system using extinguishers directly in the electronics racks that are activated using a serie of smoke detectors. In that case we will recommend powder extinguisher as they can be used in the electronics and at the same time they can not produce a problem related with air displacement like the CO_2 extinguishers.

7.5 Cryogenics

For the recovery of the Xenon to the recovery bottles the use of liquid Nitrogen is needed. The recovery bottles are in a thermally insulated container that on one hand makes the process more efficient and also prevents Nitrogen to spill in the experimental area.

When using liquid Nitrogen individual protection such as gloves and face protection is required.

8 Risk assesment for gas

The appendix describes the riks assesment of gas in the experimental area for the HARP magnet during the operation of the NEXT-DEMO detector.

8.1 Characteristics of the area

The location of the experimental area can be seen in Fig. 12. The size and distribution of the area can be seen in 13. The area has a total surface of $25m^2$, and a general view of the area is presented in Fig. 14.

On of the walls of the area is a metalic non tight wall with a gate. The rest of the walls and the floor seem relatively gas tight.

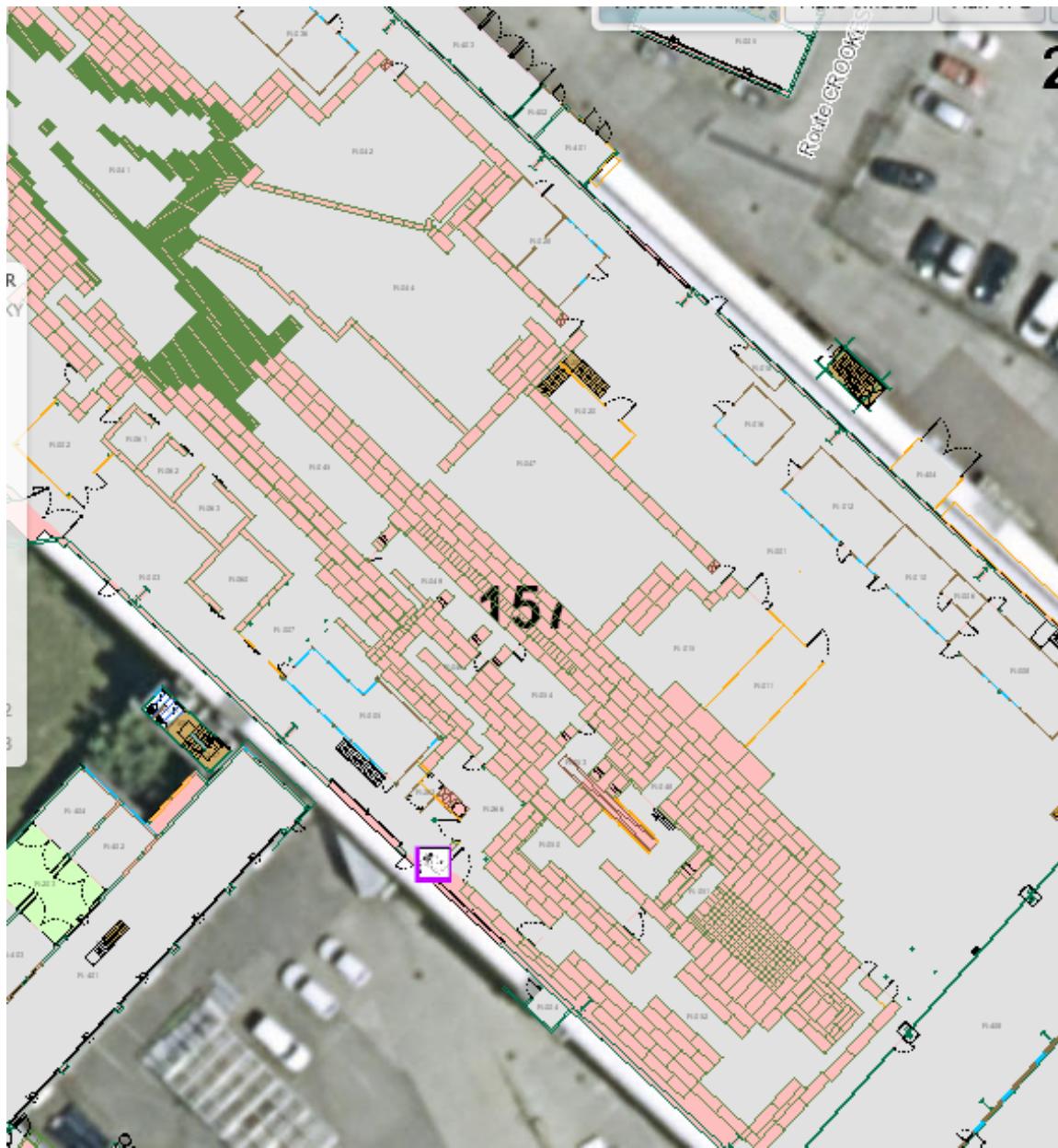


Figure 12. Location of the experimental area

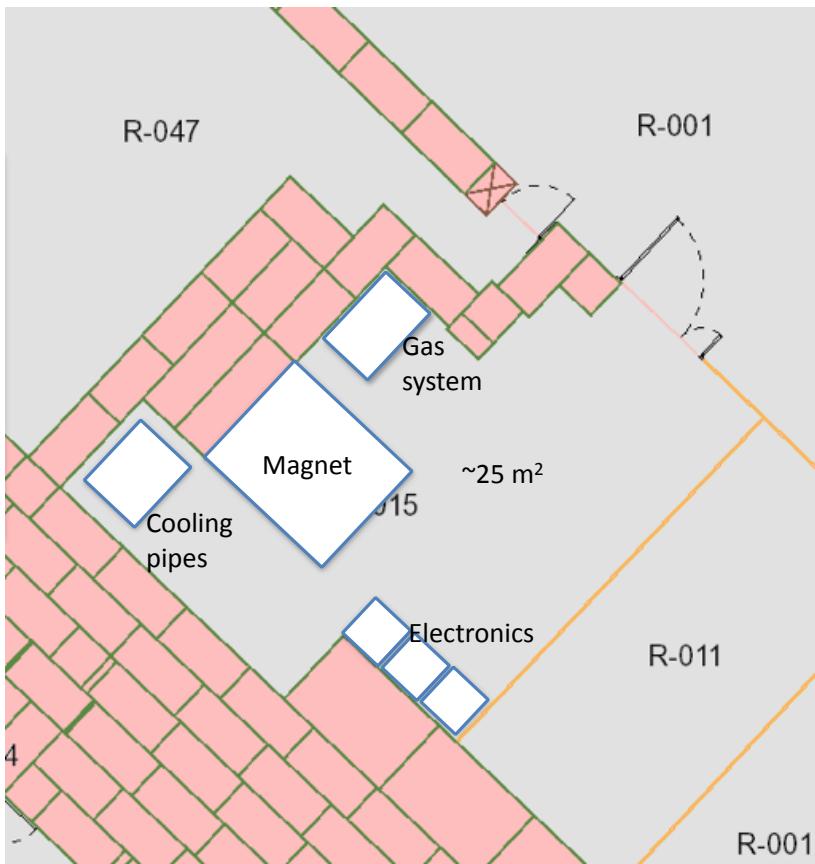


Figure 13. Proposed distribution of the NEXT-DEMO, magnet and all different necessary systems for the operation of the detector in the experimental area



Figure 14. Proposed distribution of the NEXT-DEMO, magnet and all different necessary systems for the operation of the detector in the experimental area

8.2 Related Gas systems

Only one gas system will be used in the area. The gas system has been described in the main document. The following table summarizes the quantity and distribution of the different gases to be used in NEXT-DEMO:

Gas Type	Supply description	Volume (gas at ambient P ad T)	Supply location (see figure 13)	Comments
Xenon	One bottle	255 liters	GS 1	
CO2	One bottle	200 liters (approx.)	GS 1	

Table 3. Gas supply for NEXT-DEMO

8.3 Gas release risk assesment

8.3.1 Relevant Facts

- All gases are at room temperature during normal operation.
- A gas leak of the order of 1% is easily detected.
- All used gases are transparent and odorless.

8.3.2 Risk assesment

Scenario 1

A leak due to a material breaking of a material impact that causes a hole with a size of the a 1/4" in the gas system. Total loss of the gas.

Consequences: The total volume will around 15 minutes for a complete release to the atmosphere. The total volume of Xenon is described in the **Safety** section of this document. The drop of pressure will be detected by the slow control system and it will stop the pump. People should leave the area until the pressure reaches a stable level in all different pressure gauges.

Scenario 2

A leak occurs during cryogenic recovery of the Xenon.

Consequences: In this case the slow control will be innactivated. In the other hand, most of the Xenon should be frozen and the total volume of Xenon released to the atmosphere won't be relevant.

8.4 Conclusions

8.4.1 Oxygen detection hazard

Takin in consideration all possible risks and the mitigation actions, it is recommended to install Oxygen Detection Heads.

Local flashing lights inside the area as well as at each entry can be installed and connected to the slow control if necessary.

A List of attachments

This document is accompanied by a list of technical documents that reflect all the technical specifications and test performed in the gas system, pressure vessel and recovery bottles.

These documents are the following:

- "Pressure Vessel" folder: It is a folder with the calculations, drawings and hidrostatic test performed for the construction fo the NEXT-DEMO pressure vessel.
- "Gas System" older: It contains the technical description of the gas system (file VALCI-MONT-195.pdf), the data sheets of all different components of the gas system and the pneumatic test of the system.
- "Recovery bottle" folder: It contains the information about the test of bottles for the cryo recovery.
- "Cables and connections FUG" file: It provides information about the cables used for the high voltage.
- "HAMEG-MAN-DE-EN-HMPSeries" file: It is the tecnical information for the power supply of the electronics.
- "HHV Power Supplies" folder: It contains all the technical information related with the high voltage modules.